

Steens Physical®

UNIVERSAL TREATMENT TABLE



RELIABLE, VERSATILE AND AFFORDABLE

The Universal Treatment Table from Steens Physical®:
A precision tool for chiropractors and manual therapists

This is one of our most advanced products - developed over several years in cooperation with leading therapists in Norway and Denmark. It has undergone continuous development and is now considered to be one of the best treatment tables on the market. This table is constructed to offer the user a multitude of treatment modalities through a variety of drops and adjustable functions.



Features

Steens Physical® Universal Treatment Table combines unique durability, a multitude of treatment options and a competitive price compared to therapeutic devices with similar features.

- Drop function can be built into the thoracic section and the pelvic section if required
- Length adjustable ankle support with traction belt
- Flexion with lateral flexion. Both free movement and locked position of the pelvic section in flexion and lateral position



Developed in cooperation with chiropractors and manual therapists


Our Universal Treatment Table is developed in close cooperation with the Norwegian chiropractors Knut Assjer and Vilhelm Scheel who introduced the flexion distraction concept.

Assjer and Scheel were pioneers in developing the first Flexion Distraction Universal Treatment Table in Norway, presenting a therapeutic tool that gives the unique possibility of combining all basic manual techniques.

Steens Physical® Universal Treatment Table is based on these principles and is more polifacetic than other similar tables. The patient can be treated in dorsal-, ventral- and side postures.



Electric motor control for smooth and precise height adjustments between 50 and 93 cms.

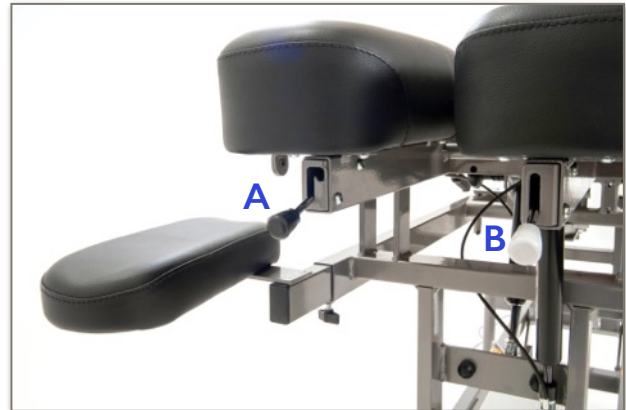

A: Control for thoracic section release
B: Control for gas suspended adjustable head section



Extra drops available as option



Flexion handle and ankle belt



Adjustable head and arm supports

Distributor

Technical specifications:

Length: 174 - 190 cm
Width: 56 cm
Max height: 93 cm
Min height: 50 cm
Weight: 65 kg
Lifting capacity: 130 kg
Caring capacity: 250 kg